

ErgoTable



- Ergonomic working
- Sophisticated design
- Sturdy construction



Ergologically!

Intelligent construction for greater efficiency.



The CSSD packing station is a complex working space with high demands on the employee. The ErgoTable height adjustable packing table has been developed by MMM to alleviate the sometimes physically demanding work and meets all the requirements of ergonomics, hygiene and modular design. The working height is adjustable to the required level from 800 mm to 1050 mm by pressing a button; the back is protected regardless of whether working in a sitting or standing position. Thanks to the ergonomic circular cut-out at the front of the tabletop the high mounted bins and shelves are easy and straightforward to access.

The well illuminated workspace enables effective visual and functional testing of technically complex medical devices and with the many useful compo-

nents such as hook racks, wire mesh bins and cutting devices, all the necessary utensils for screen and instrument packing are within easy reach.

The ErgoTable is kept particularly tidy and hygienic through the media duct, which runs under the rear tabletop. Digital and power cables are stored here for connecting a monitor, keyboard and other equipment and fed to the central cable outlet. The monitor mount is adjustable horizontally and vertically to the optimal viewing angle and can even hold a mini PC. The clean design solution and the polished stainless steel finish of the modules make the ErgoTable both completely hygienic as well as ergonomic.

Technical values

Dimensions mm	Working height mm	Material	Colours	Electrical connection
2060 x 900	800 - 1050	HPL	Capri blue	Schuko socket 230 V/50 Hz AC
2060 x 900	800 - 1050	HPL	Off-white	Schuko socket 230 V/50 Hz AC

Other colours on request

Ergonomic

Many details make the ErgoTable a true ergonomic table: Height adjustable, round cut-out, monitor position, illumination etc.

Hygienic

No more loose cabling over the packing table. Digital and power cables are tidily stored away in the pull-out media duct.

Individual

Every CSSD is different, as are the demands of the respective work routines. With the many different modules to equip the ErgoTable your individual needs are addressed.

Horizontally and vertically adjustable monitor stand.




Media duct to accommodate digital and power cables.



Simple height adjustment at the press of a button.



ErgoTable can be put together without gaps and therefore without the risk of pinching between the backs.

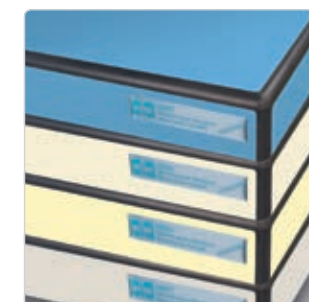

Our standard colours are Capri Blue and Off-white. Other colour options are available upon request.



Drawer pedestal with individually configurable interior.



MMM Group

MMM has been operating worldwide as one of the leading system providers in the service of health since 1954. With a complete range of products and services relating to all aspects of cleaning, disinfection and sterilization systems for the areas of **Healthcare** and **Life Science**, MMM has positioned itself as a crucial quality and innovation driver in the German and international

market. Our products are individually adapted to the requirements of our customers all over the world. The high vertical range of manufacturing in our production plants ensures that we fulfil the strictest demands of quality in the medical technology sector. More than 1100 employees apply their expertise and dedication to the mission of the MMM Group: **Protecting human health.**



Health
made in Germany

MMM Münchener
Medizin Mechanik GmbH
Simmelweisstraße 6
82152 Planegg / Munich, Germany
Phone: +49 89 89918-0
Fax: +49 89 89918-118
www.mmmgroup.com

